

COMPASS Synoptic CB Porewater: DOC

July 2023 Samples

2025-10-31

```
##Setup - Change things here & write any notes
```

```
#identify section  
cat("Setup Information")
```

```
## Setup Information
```

```
##### Run information - PLEASE CHANGE  
Date_Run = "07/31/23" #Date that instrument was run  
Run_by = "Stephanie J. Wilson" #Instrument user  
Script_run_by = "Stephanie J. Wilson" #Code user  
run_notes = "NA" #any notes from the run  
samples <- c("GCW", "GWI", "MSM", "SWH") #whatever identifies your samples within the same names  
samples_pattern <- paste(samples, collapse = "|")  
#samples_pattern <- "GCW" #use this instead of the line above if you have only one site code  
chks_name = "Chk_Std_50ppmC_2ppmN" #what did you name your check standards?
```

```
##### File Names - PLEASE CHANGE
```

```
#file path and name for raw summary data file  
#raw_file_name = "Raw Data/COMPASS_SynopticCB_PW_DOC_202505.txt" #example  
raw_file_name = "Raw Data/TOCTN_COMPASS_Synoptic_TCTN_202307.txt"
```

```
#file path and name for raw all peaks file
```

```
#raw_allpeaks_name = "Raw Data/COMPASS_SynopticCB_PW_DOC_202505_allpeaks.txt" #example  
# raw_allpeaks_name = "Raw Data/COMPASS_SynopticCB_PW_DOC_202411_allpeaks.txt"
```

```
#file path and name of processed data file
```

```
#processed_file_name = "Processed Data/COMPASS_SynopticCB_PW_Processed_DOC_202505.csv" #example  
processed_file_name = "Processed Data/COMPASS_SynopticCB_PW_Processed_DOC_202307.csv" #example
```

```
##### Log Files - PLEASE CHECK
```

```
#downloaded metadata csv - downloaded from Google drive as csv for this year  
Raw_Metadata = "Raw Data/COMPASS_SynopticCB_PW_SampleLog_2023.csv"
```

```
#qaqc log file path for this year
```

```
#Log_path = "Raw Data/COMPASS_Synoptic_TOCTN_QAQClog_2024.csv"
```

```
##Set Up Code
```

```
##Read in metadata and create similar sample IDs for matching to samples
```

Import Data Functions

Import Sample Data

```
## Import Sample Data

## New names:
## * ` ` -> '...14'

## # A tibble: 6 x 4
##   sample_name      npoc_raw tdn_raw run_datetime
##   <chr>          <dbl>    <dbl>    <chr>
## 1 GCW_202307_UP_LysA_20cm 49.1     2.08    7/31/2023 4:45:43 AM
## 2 GCW_202307_TR_LysA_20cm 14.0     0.371   7/31/2023 5:07:12 AM
## 3 GCW_202307_TR_LysA_45cm 44.3     2.10    7/31/2023 5:36:14 AM
## 4 GCW_202307_TR_LysB_10cm 11.4     0.410   7/31/2023 6:01:31 AM
## 5 GCW_202307_TR_LysB_45cm 4.38     0.458   7/31/2023 6:23:05 AM
## 6 GCW_202307_TR_LysC_10cm -0.0478  0.0385  7/31/2023 10:49:23 AM
```

Assessing Standard Curves - assessed manually on the instrument


Assess Check Standards

```
## Assess the Check Standards

## New names:
## * ` ` -> '...14'

## [1] "Carbon CHECK STANDARD RSD TOO HIGH - REASSESS"

## [1] "Nitrogen CHECK STANDARD RSD TOO HIGH - REASSESS"
```

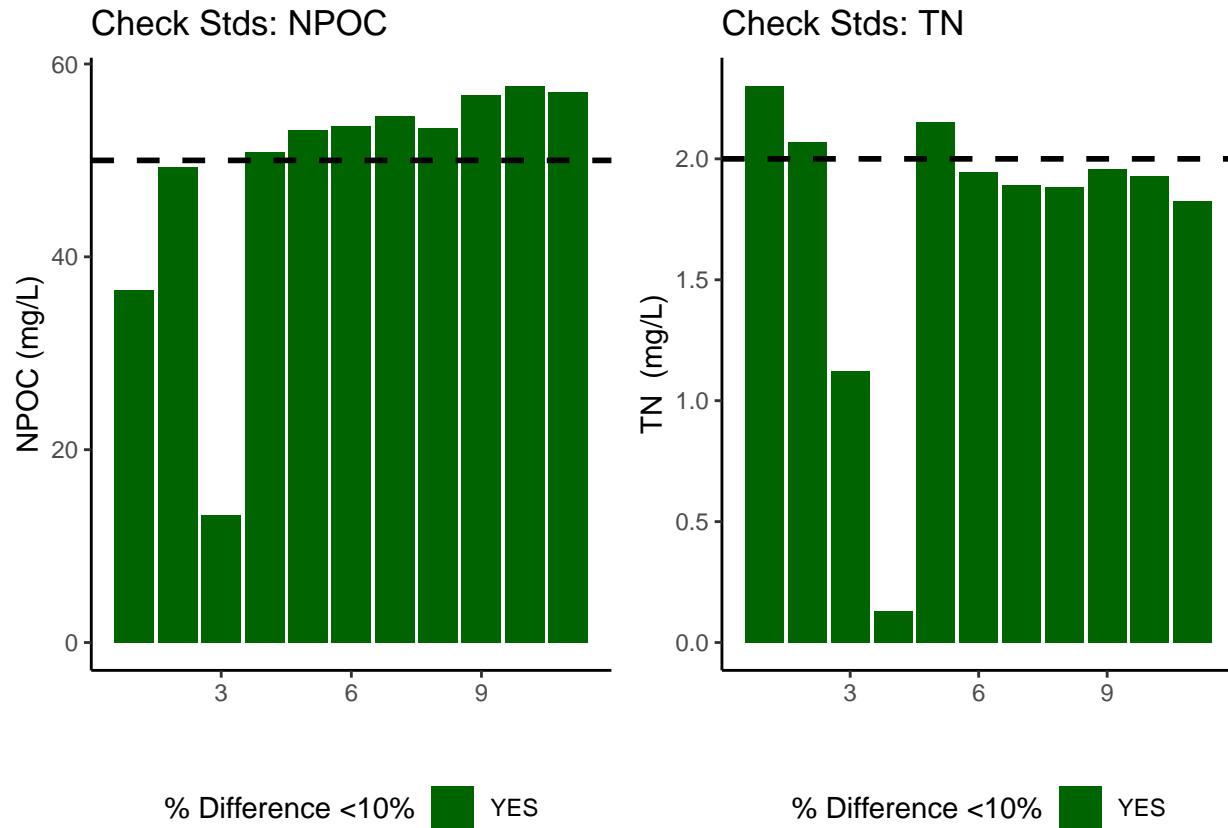


```
## [1] ">60% of Carbon Check Standards are within range of expected concentration"
## [1] ">60% of Nitrogen Check Standards are within range of expected concentration"
```

Assess Blanks

```
## Assess Blanks

## New names:
## * `` -> '...14'

## [1] ">60% of Carbon Blank concentrations are lower 25% quartile of samples"
## [1] ">60% of Nitrogen Blank concentrations are lower 25% quartile of samples"
```

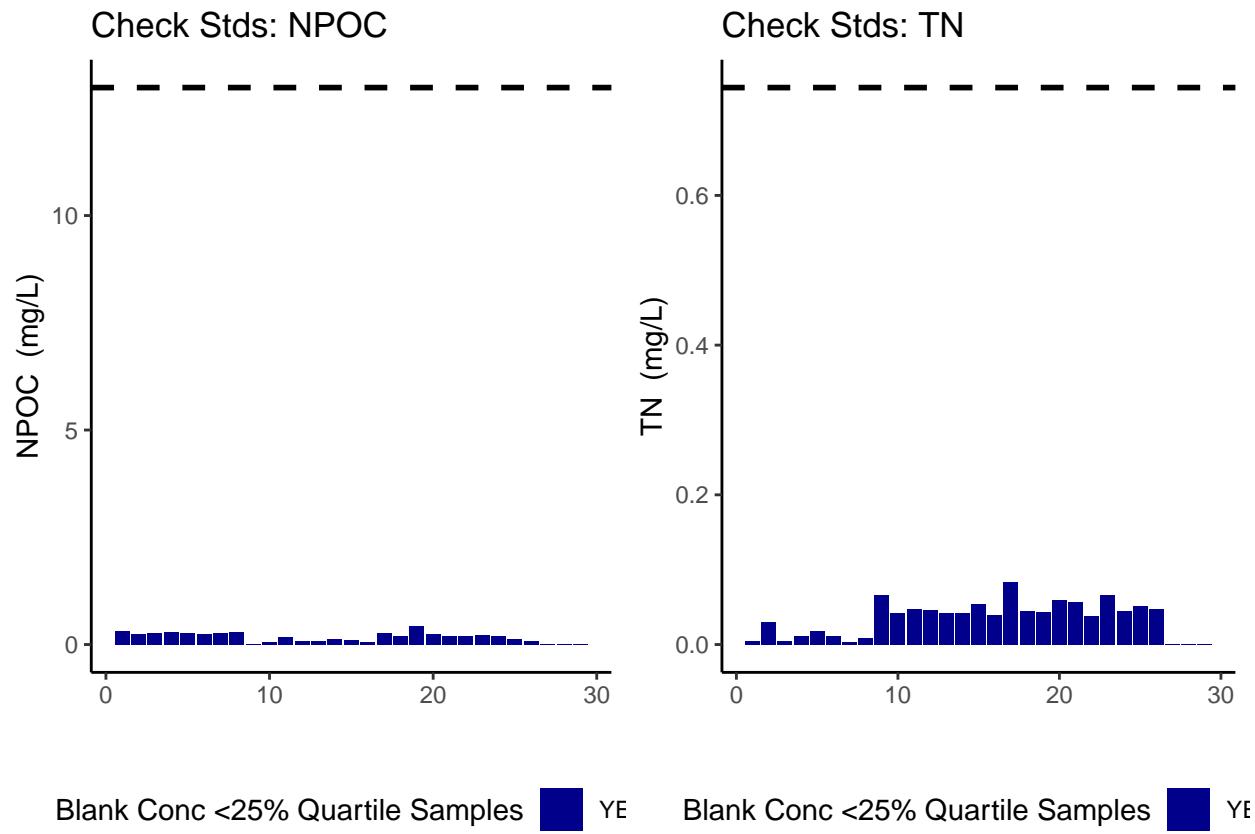

```
## carbon blanks:
```

```
## [1] 0.1724955
```


```
## nitrogen blanks:
```

```
## [1] 0.03486172
```

Assess Duplicates

```
## Assess Duplicates
```

```
## Warning: Using `size` aesthetic for lines was deprecated in ggplot2 3.4.0.
## i Please use `linewidth` instead.
## This warning is displayed once every 8 hours.
## Call `lifecycle::last_lifecycle_warnings()` to see where this warning was
## generated.
```

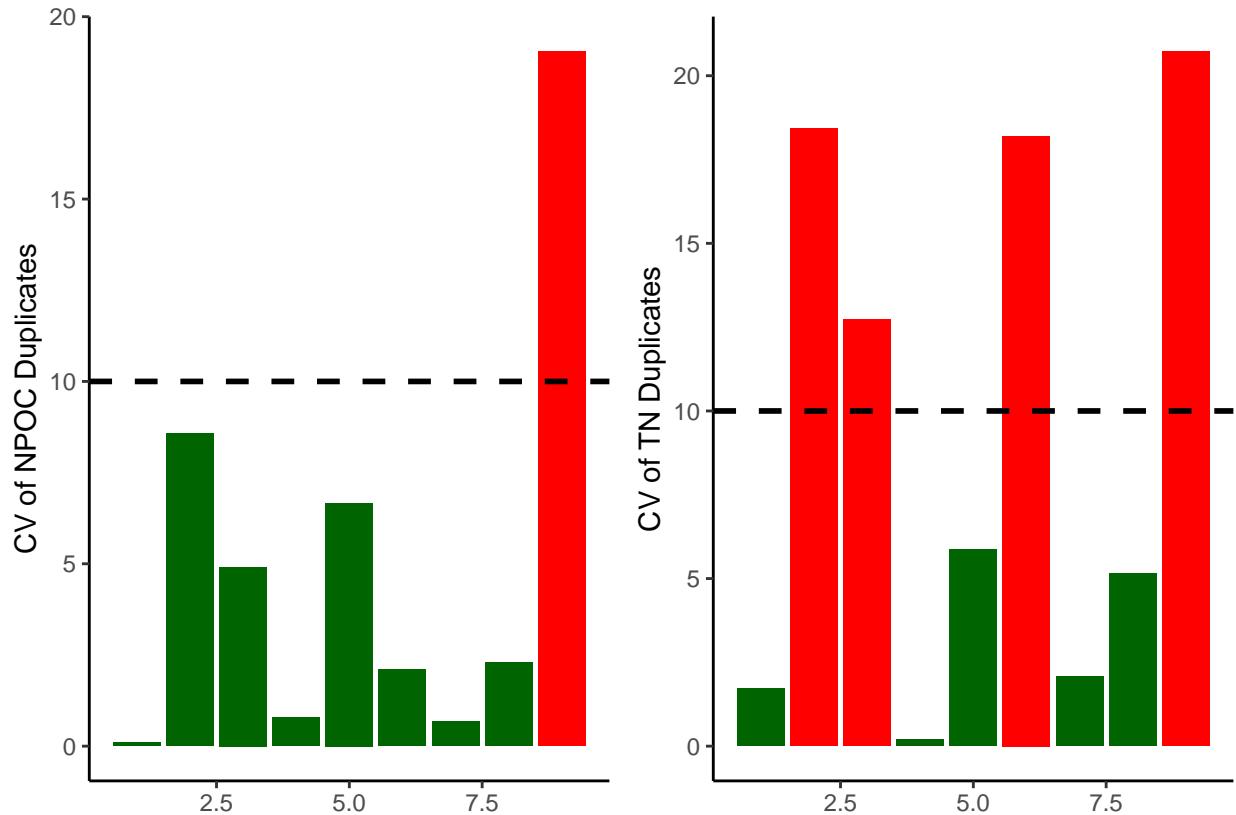


```
## [1] ">60% of Carbon Duplicates have a CV <10%"
```

```
## [1] ">60% of Nitrogen Duplicates have a CV <10%"
```


Sample Flagging - Are samples Within the range of the curve?

```
## Sample Flagging
```


C: Grey = Within Range of Curve

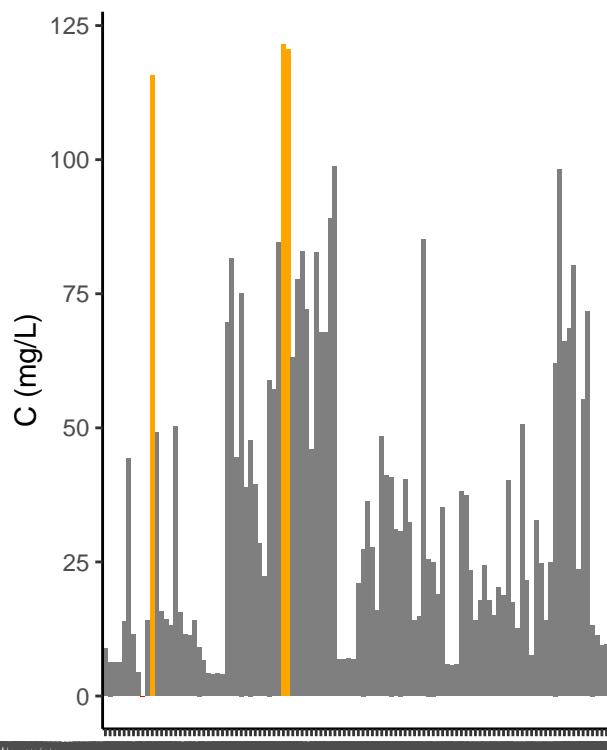

Samples: NPOC

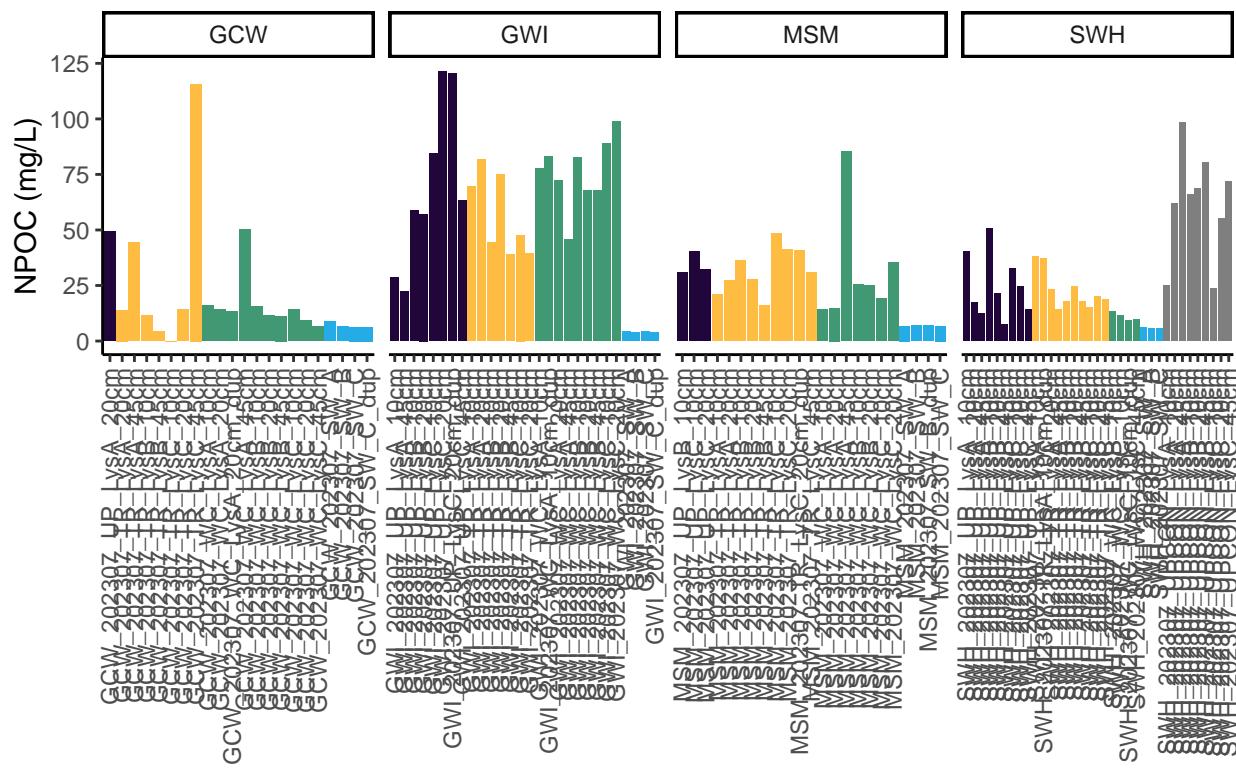

Samples: TN

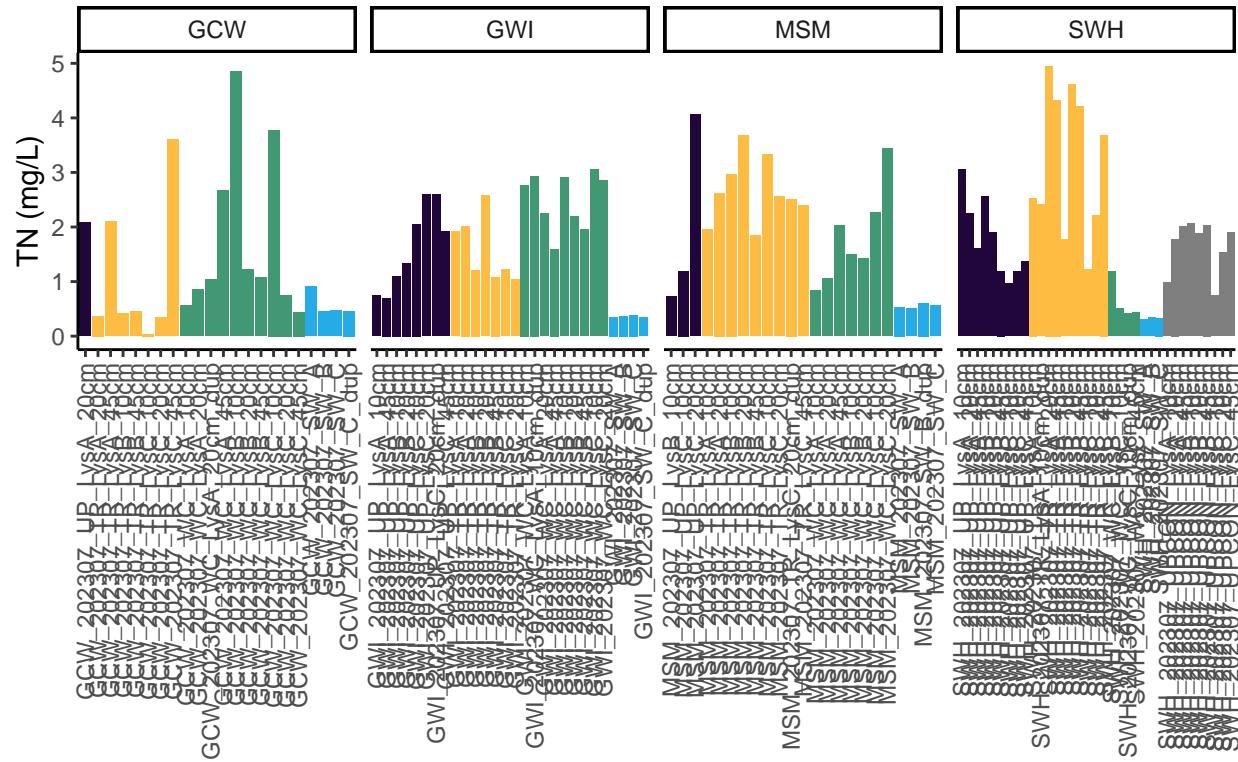

Convert data from mg/L to uMoles/L

Check to see if samples run match metadata & merge info

```
## Check Sample IDs with Metadata
## All sample IDs are present in metadata.
```

Export Processed Data

```
## Export Processed Data

## # A tibble: 6 x 21
##   Project     Region Site Zone Replicate Depth_cm Sample_ID Year Month Day
##   <chr>      <chr>  <chr> <fct> <chr>       <int> <chr>    <int> <int> <int>
## 1 COMPASS: Sy~ CB    GCW   UP    A           20  GCW_2023~ 2023    7    11
## 2 COMPASS: Sy~ CB    MSM   UP    B           10  MSM_2023~ 2023    7    13
## 3 COMPASS: Sy~ CB    MSM   UP    C           10  MSM_2023~ 2023    7    13
## 4 COMPASS: Sy~ CB    MSM   UP    C           20  MSM_2023~ 2023    7    13
## 5 COMPASS: Sy~ CB    SWH   UP    A           10  SWH_2023~ 2023    7    19
## 6 COMPASS: Sy~ CB    SWH   UP    A           20  SWH_2023~ 2023    7    19
## # i 11 more variables: Time <lgl>, Time_Zone <lgl>, npoc_mgL <dbl>,
## #   npoc_uM <dbl>, npoc_flag <chr>, tdn_mgL <dbl>, tdn_uM <dbl>,
## #   tdn_flag <chr>, Analysis_runtime <chr>, Run_notes <chr>, Field_notes <chr>
```

#end