

# Lucie Stetten Incubation Samples: SO<sub>4</sub>/Cl

2025 Samples

2025-12-01

## Contents

0.1	Run Information . . . . .	2
0.2	Assess Standard Curves . . . . .	3
0.3	Assess Check Standards . . . . .	4
0.4	Assess Blanks . . . . .	5
0.5	Assess Duplicates . . . . .	7
0.6	Calculate mmol/L concentrations & salinity, add dilutions . . . . .	8
0.7	Assess Analytical Spikes . . . . .	9
0.8	Check if samples within the range of the standard curve . . . . .	10
0.9	Export Processed Data . . . . .	10

##Add Required Packages

## 0.1 Run Information

```
##### Run information - PLEASE CHANGE
Date_Run = "2025-11-24" #Date that instrument was run
Run_by = "Zoe Read" #Instrument user
Script_run_by = "Zoe Read" #Code user
run_notes = "No issues.
" #any notes from the run
samples <- c("GCW", "GWI", "MSM", "SWH") #whatever identifies your samples within the same names
samples_pattern <- paste(samples, collapse = "|")
#samples_pattern <- "GCW" #use this instead of the line above if you have only one site code

##Fix samples that were entered incorrectly
samples_to_remove <- NA
Old_ID_1 = NA
New_ID_1 = NA

##### File Names - PLEASE CHANGE
#file path and name for raw summary data file
raw_file_name_cl = "Raw Data/COMPASS_Lucies_Samples_20251124_Cl.txt"
raw_file_name_so4 = "Raw Data/COMPASS_Lucies_Samples_20251124_S04.txt"

#file path and name of processed data file
processed_file_name = "Processed Data/COMPASS_Lucies_Samples_Processed_Cl_S04_20251124.csv"

##### Log Files - PLEASE CHECK
#downloaded metadata csv - downloaded from Google drive as csv for this year
Raw_Metadata = "Raw Data/COMPASS_SynopticCB_PW_SampleLog_2025.csv"

#qaqc log file path for this year
Log_path = "Raw Data/COMPASS_Synoptic_Cl_S04_QAQClog_2024.csv"
```


##Set Up Code - constants and QAQC cutoffs

##Import Sample Data


##Fix samples entered incorrectly

## 0.2 Assess Standard Curves


Chloride Std Curve

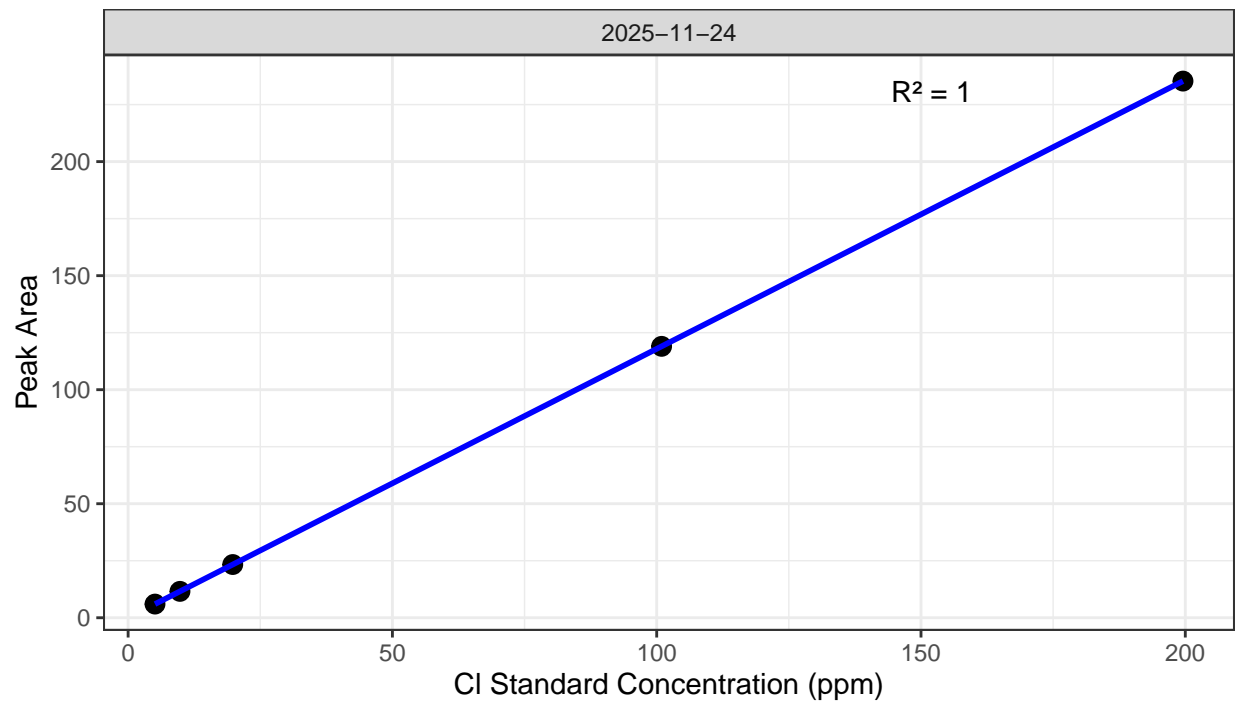

Sulfate Std Curve

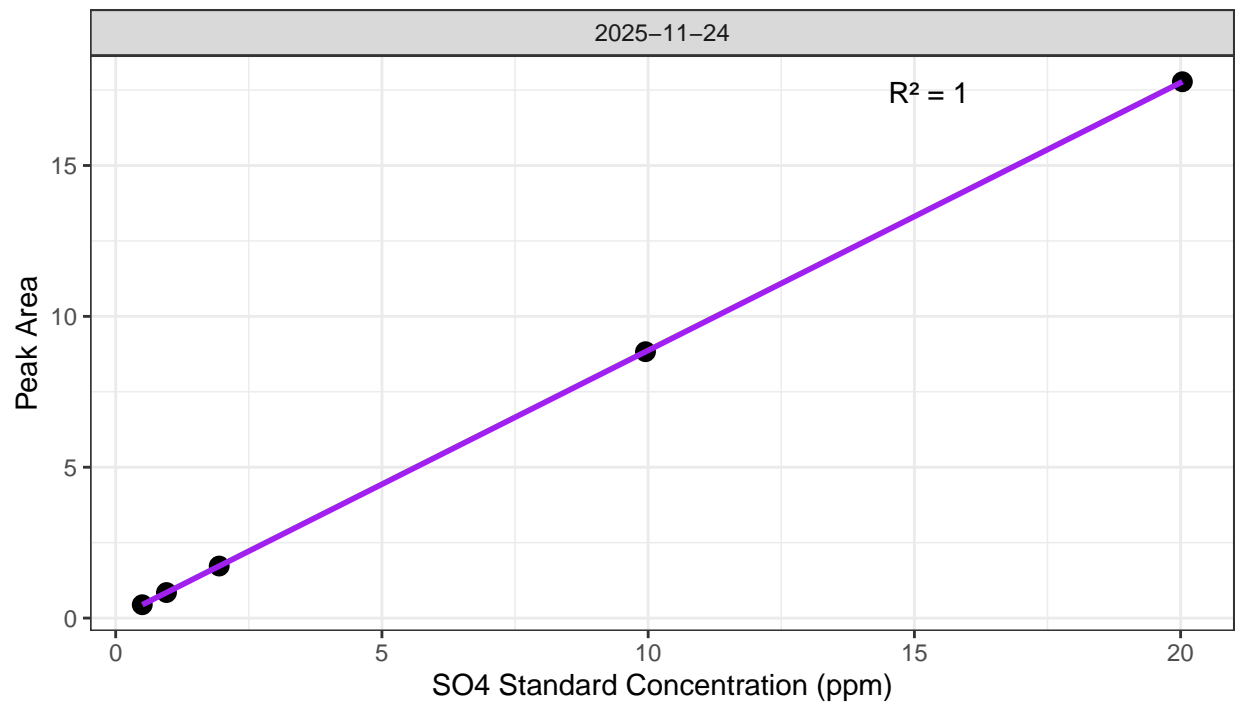

```
## [1] "QAQC log file exists and has been read into the code."
```

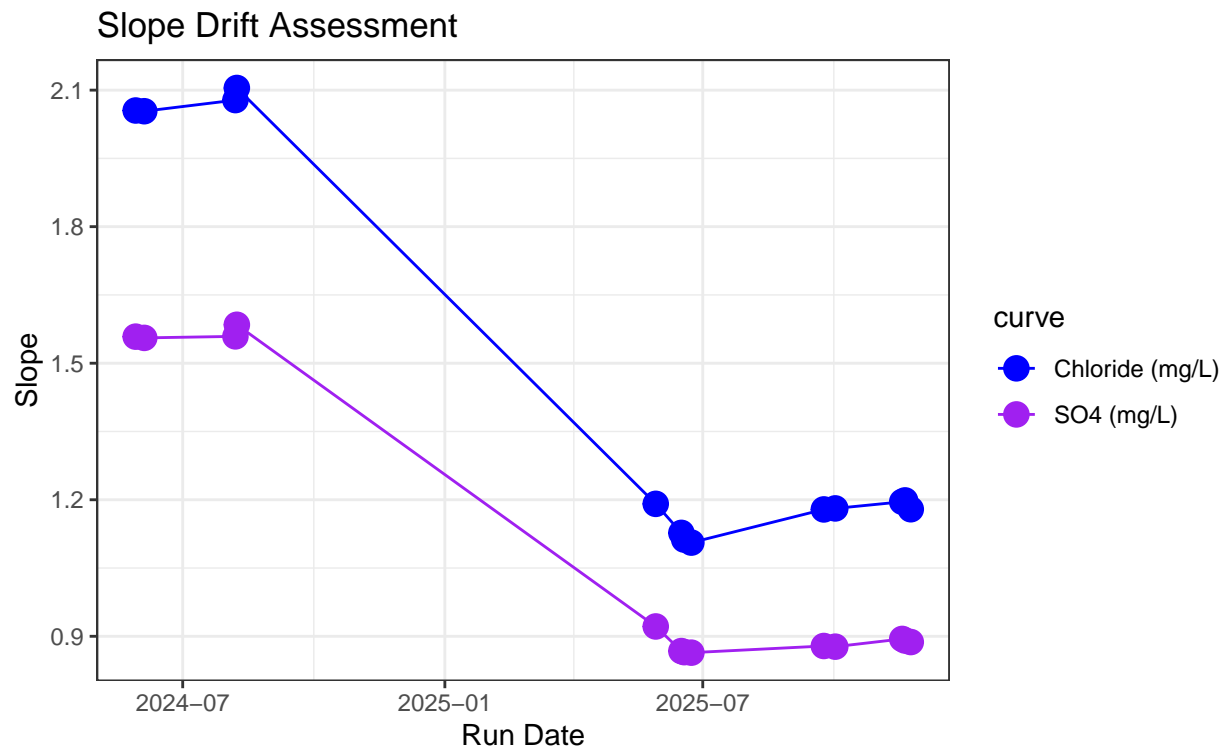

```
## [1] "Cl Curve r2 GOOD"
```

```
## [1] "SO4 Curve r2 GOOD"
```


### 0.3 Assess Check Standards

```
## # A tibble: 5 x 5
##   sample_ID mean_Cl sd_Cl cv_Cl flag_Cl
##   <chr>      <dbl> <dbl> <dbl> <chr>
## 1 Standard 1 5.19 0.0702 0.0135 Chloride Check Standard RSD within Range --
## 2 Standard 2 9.99 0.120 0.0120 Chloride Check Standard RSD within Range --
## 3 Standard 3 20.0 0.149 0.00744 Chloride Check Standard RSD within Range --
## 4 Standard 4 101. 0.368 0.00363 Chloride Check Standard RSD within Range --
## 5 Standard 5 199. 0.141 0.000708 Chloride Check Standard RSD within Range --
```

```
## # A tibble: 5 x 5
##   sample_ID mean_SO4 sd_SO4 cv_SO4 flag_SO4
##   <chr>      <dbl> <dbl> <dbl> <chr>
## 1 Standard 1 0.504 0.00307 0.00610 Sulfate Check Standard RSD within Range --
## 2 Standard 2 0.959 0.00345 0.00359 Sulfate Check Standard RSD within Range --
## 3 Standard 3 1.95 0.00697 0.00357 Sulfate Check Standard RSD within Range --
## 4 Standard 4 10.0 0.0478 0.00478 Sulfate Check Standard RSD within Range --
## 5 Standard 5 20.1 0.0219 0.00109 Sulfate Check Standard RSD within Range --
```

```
## [1] ">80% of Chloride Check Standards have RSD within range - PROCEED"
```

```
## [1] ">80% of Sulfate Check Standards have RSD within range - PROCEED"
```

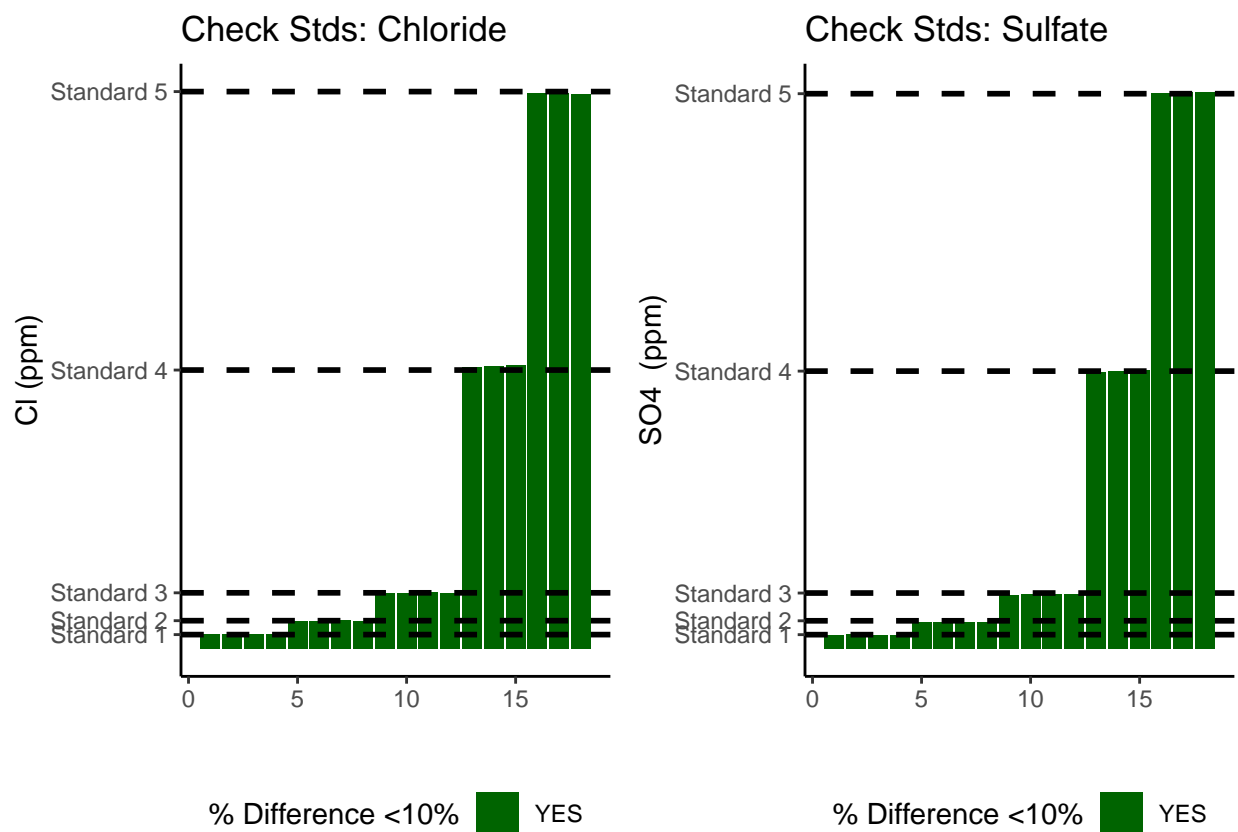


```
## [1] ">80% of Chloride Check Standards are within range of expected concentration - PROCEED"
```

```
## [1] ">80% of Sulfate Check Standards are within range of expected concentration - PROCEED"
```

#### 0.4 Assess Blanks

```
## [1] ">80% of Chloride Blank concentrations are lower 25% quartile of samples"
```

```
## [1] ">80% of Sulfate Blank concentrations are lower 25% quartile of samples"
```

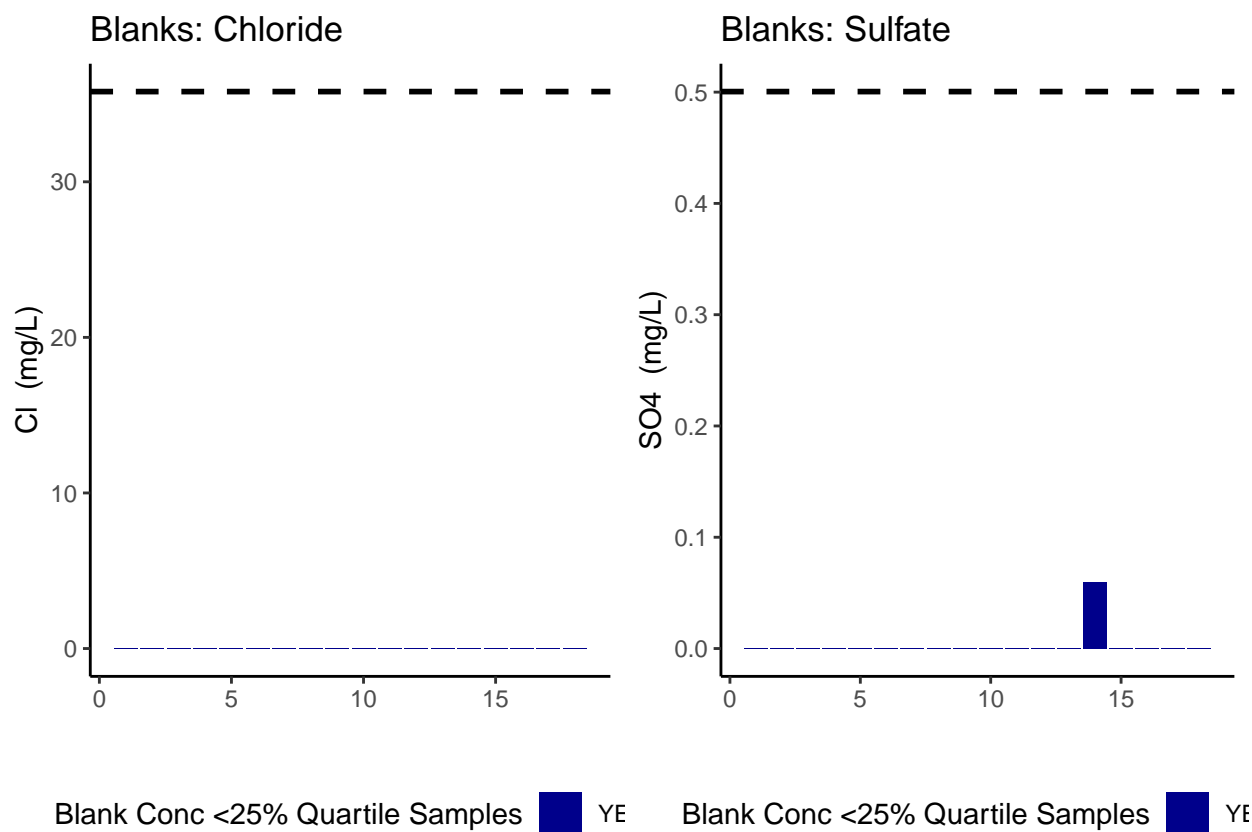


```
## Chloride blanks mean ppm:
```

```
## [1] 0.006422222
```

```
## Sulfate blanks mean ppm:
```


```
## [1] 0.003327778
```

## 0.5 Assess Duplicates

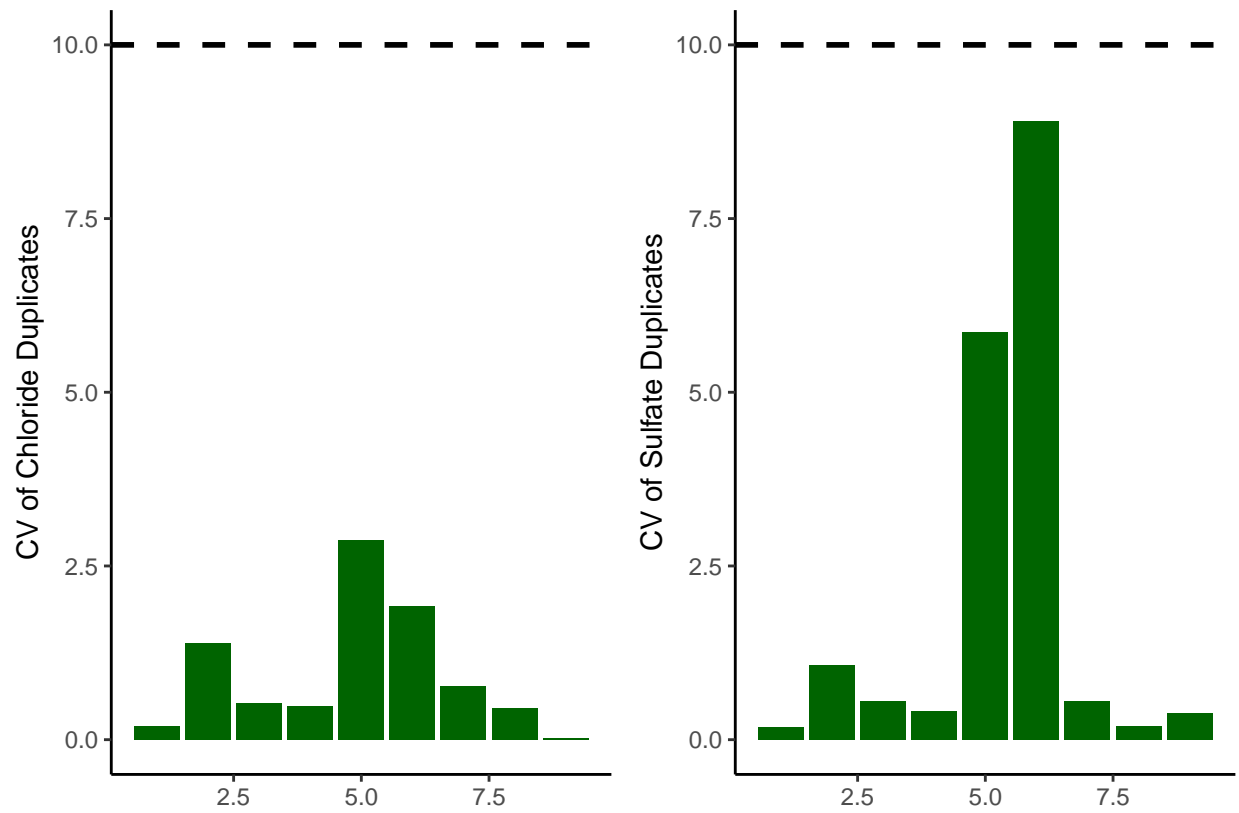

```
## [1] ">80% of Chloride Duplicates have a CV <10% - PROCEED"
```

```
## [1] ">80% of Sulfate Duplicates have a CV <10% - PROCEED"
```

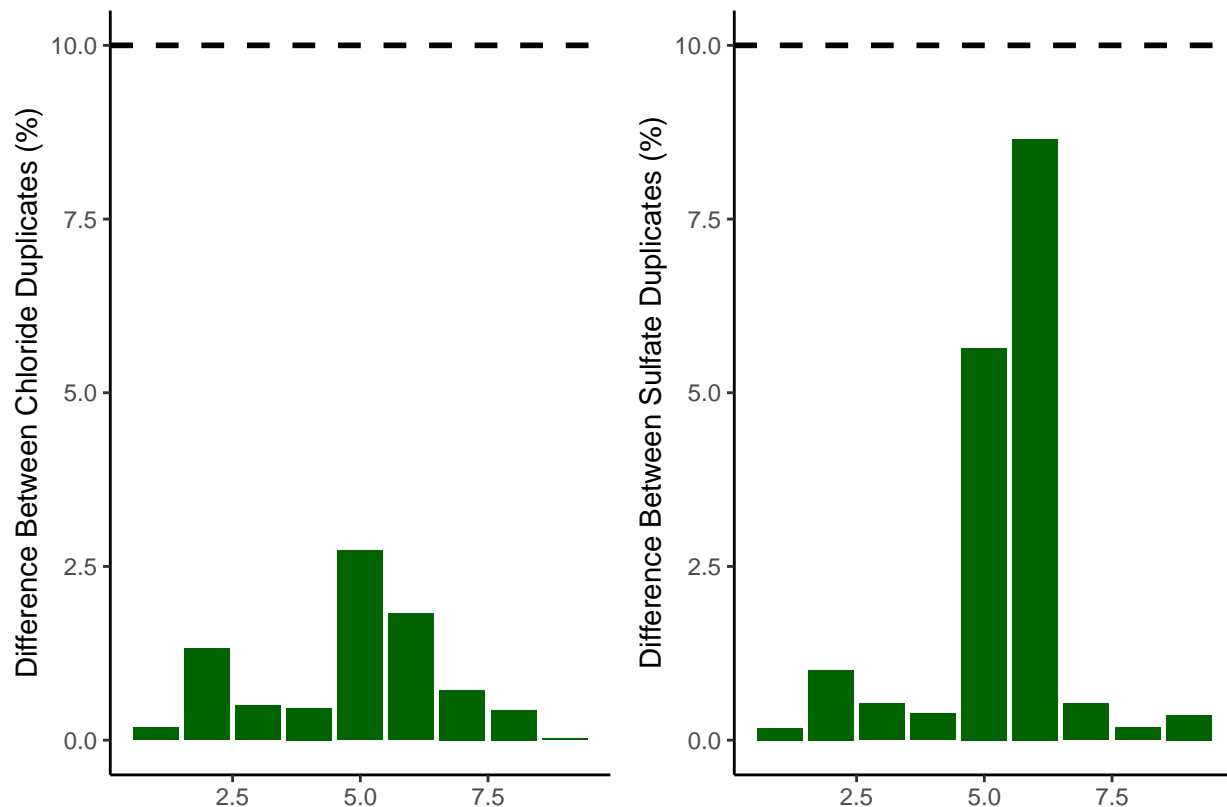

```
## [1] ">80% of Chloride Duplicates have a percent difference <10% - PROCEED"
```

```
## [1] ">80% of Sulfate Duplicates have a percent difference <10% - PROCEED"
```

## 0.6 Calculate mmol/L concentrations & salinity, add dilutions

```
# Convert ppm to mmol/L
all_dat$SO4_Conc_mM <- (all_dat$SO4_ppm / s_mw)
all_dat$Cl_Conc_mM <- (all_dat$Cl_ppm / cl_mw)

# Calculate Salinity
# calculated using the Knudsen equation
# Salinity = 0.03 + 1.8050 * Chlorinity
# Ref: A Practical Handbook of Seawater Analysis by Strickland & Parsons (P. 11)
# =((1.807*Cl_ppm)+0.026)/1000
all_dat$salinity <- ((1.8070 * all_dat$Cl_ppm) + 0.026) / 1000

#Need to determine dilution factors for your samples
#for Steph / COMPASS this depends on the site so...
all_dat$Dilution <- 1
all_dat$Dilution <- ifelse(str_detect(all_dat$sample_ID, "MSM") & str_detect(all_dat$sample_ID, "UP"),
all_dat$Dilution <- ifelse(str_detect(all_dat$sample_ID, "MSM") & str_detect(all_dat$sample_ID, "TR"),
all_dat$Dilution <- ifelse(str_detect(all_dat$sample_ID, "MSM") & str_detect(all_dat$sample_ID, "WC"),
```



```

all_dat$Dilution <- ifelse(str_detect(all_dat$sample_ID, "MSM") & str_detect(all_dat$sample_ID, "SW"),
all_dat$Dilution <- ifelse(str_detect(all_dat$sample_ID, "GCW") & str_detect(all_dat$sample_ID, "UP"),
all_dat$Dilution <- ifelse(str_detect(all_dat$sample_ID, "GCW") & str_detect(all_dat$sample_ID, "TR"),
all_dat$Dilution <- ifelse(str_detect(all_dat$sample_ID, "GCW") & str_detect(all_dat$sample_ID, "WC"),
all_dat$Dilution <- ifelse(str_detect(all_dat$sample_ID, "GCW") & str_detect(all_dat$sample_ID, "SW"),

all_dat$Dilution <- ifelse(str_detect(all_dat$sample_ID, "GWI") & str_detect(all_dat$sample_ID, "UP"),
all_dat$Dilution <- ifelse(str_detect(all_dat$sample_ID, "GWI") & str_detect(all_dat$sample_ID, "TR"),
all_dat$Dilution <- ifelse(str_detect(all_dat$sample_ID, "GWI") & str_detect(all_dat$sample_ID, "WC"),
all_dat$Dilution <- ifelse(str_detect(all_dat$sample_ID, "GWI") & str_detect(all_dat$sample_ID, "SW"),


all_dat$Dilution <- ifelse(str_detect(all_dat$sample_ID, "SWH"), 50, all_dat$Dilution)

all_dat$Dilution <- ifelse(str_detect(all_dat$sample_ID, "Lucie"), 5, all_dat$Dilution)

# head(all_dat)

```

## 0.7 Assess Analytical Spikes

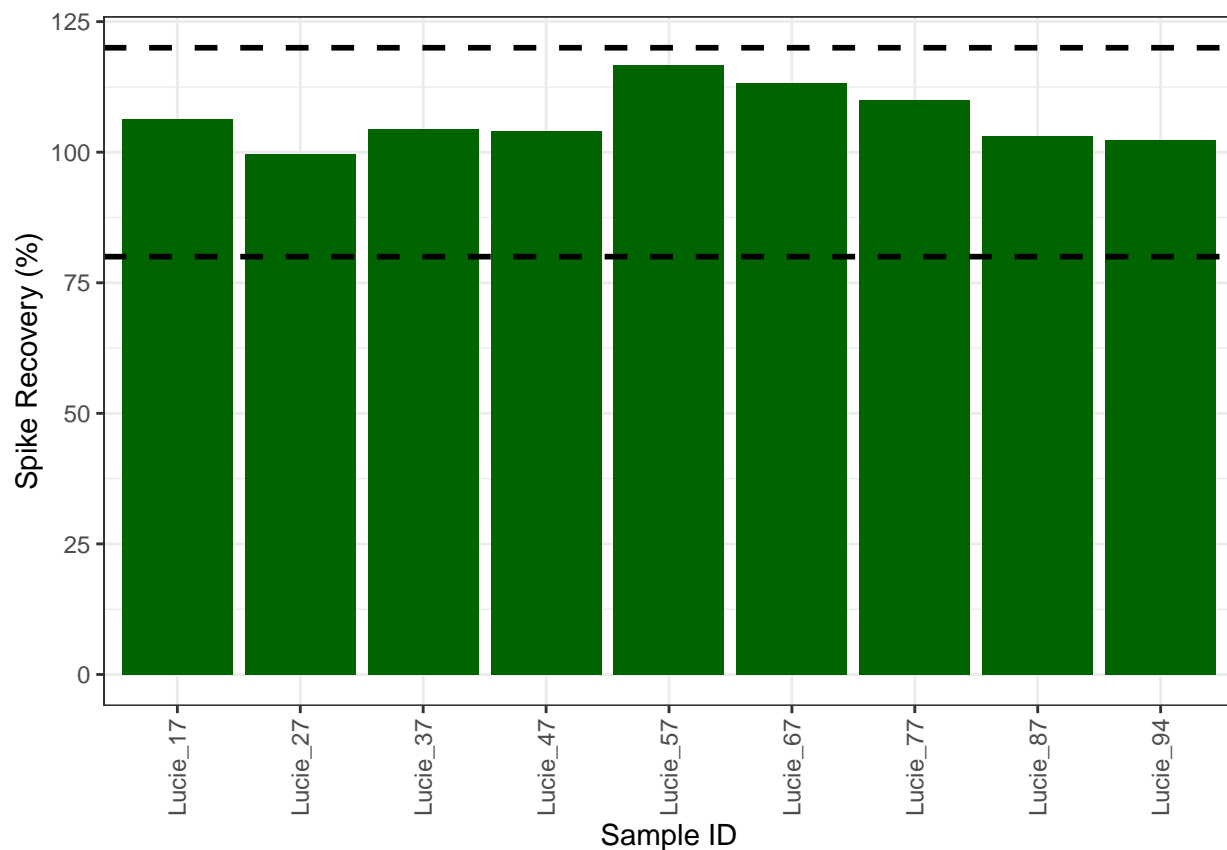

```
## [1] ">80% of S04 spikes have a recovery between the high and low cutoff - PROCEED"
```

## 0.8 Check if samples within the range of the standard curve

## Sample Flagging

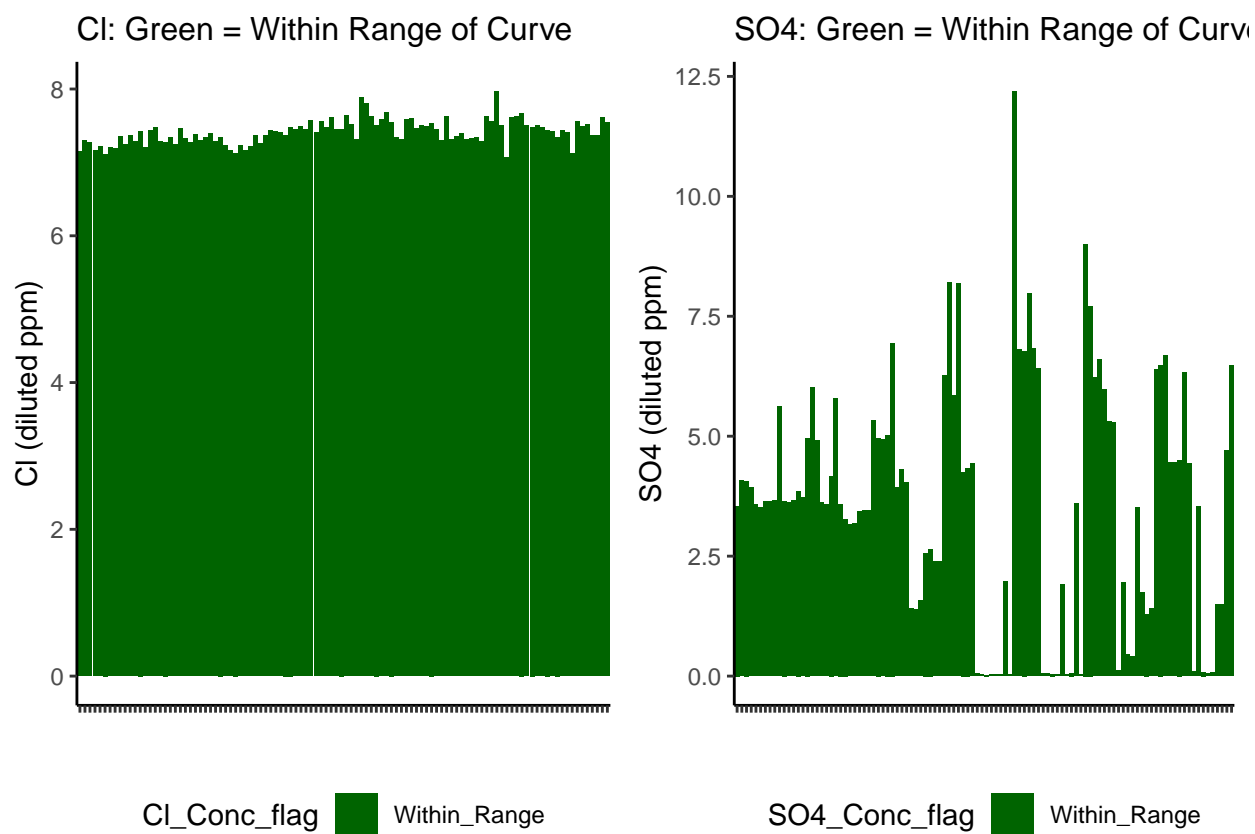


Table 1: SO4 samples

SO4_Conc_flag	Percent_samples
Within_Range	100

Table 2: Cl samples

Cl_Conc_flag	Percent_samples
Within_Range	100

## 0.9 Export Processed Data

#end