

Synoptic CB: Porewater Nutrients

November 2024 Samples

2025-08-11

Contents

0.1	Import Data & Clean	3
0.2	Assessing standard Curves	3
0.3	Dilution Corrections - ensure the latest dilution is kept	7
0.4	Performance Check	7
0.5	Analyze the Check Standards	9
0.6	Analyze Blanks	10
0.7	Analyze Duplicates	11
0.8	Spikes	12
0.9	Matrix Effects	13
0.10	Unit Converted Data Column Added (mg/L to uM)	13
0.11	Sample Flagging - Within range of standard curve	13
0.12	Pull out sample id information	13
0.13	Check to see if samples run match metadata & merge info	13
0.14	Visualize Data	14
0.15	Export Processed Data	16

##Run Information

```
cat("Run Information: NAME ") #lets you know what section you're in
```

Run Information: NAME

```
#set the run date & user name
run_date <- "08/07/2025"
sample_year <- "2024"
sample_month <- "NOVEMBER"
user <- "Isabelle Van Benschoten"

#identify the files you want to read in
#read in as a list to accommodate multiple runs in a month
NOx_files <- c("Raw Data/COMPASS_Synoptic_CB_202411_VNOx_1.csv",
              "Raw Data/COMPASS_Synoptic_CB_202411_VNOx_2.csv",
              "Raw Data/COMPASS_Synoptic_CB_202411_VNOx_3.csv")
NH3_PO4_files <- c("Raw Data/COMPASS_Synoptic_CB_202411_NH3_PO4_1.csv",
                  "Raw Data/COMPASS_Synoptic_CB_202411_NH3_PO4_2.csv",
                  "Raw Data/COMPASS_Synoptic_CB_202411_NH3_PO4_3.csv")

# Define the file path for QAQC log file - NO Need to change just check year
file_path <- "Raw Data/SEAL_COMPASS_Synoptic_QAQC_Log_2024.csv"
final_path <- "Processed Data/COMPASS_Synoptic_Nutrients_202411.csv"

#record any notes about the run or anything other info here:
run_notes <- "One PO4 check standard came back very high, all others were within range.
             A few duplicates were out of range but still within CV.
             SWH_202411_SWAMP_LysC_45cm, SWH_202411_TR_LysC_10cm,
             MSM_202411_TR_LysC_10cm, MSM_202411_TR_LysC_45cm were missing from metadata. "

#Set up file path for metadata
#downloaded metadata csv - downloaded from Google drive as csv for this year
Raw_Metadata = "Raw Data/COMPASS_SynopticCB_PW_SampleLog_2024.csv"

cat(run_notes)
```

```
## One PO4 check standard came back very high, all others were within range.
##           A few duplicates were out of range but still within CV.
##           SWH_202411_SWAMP_LysC_45cm, SWH_202411_TR_LysC_10cm,
##           MSM_202411_TR_LysC_10cm, MSM_202411_TR_LysC_45cm were missing from metadata.
```

##Setup

##Read in metadata and create similar sample IDs for matching to samples

0.1 Import Data & Clean


0.2 Assessing standard Curves

```
#Pull out standards data
```

```
## Assess Standard Curves
```


```
#Plot standards data
```

```
## Assess Standard Curves
```

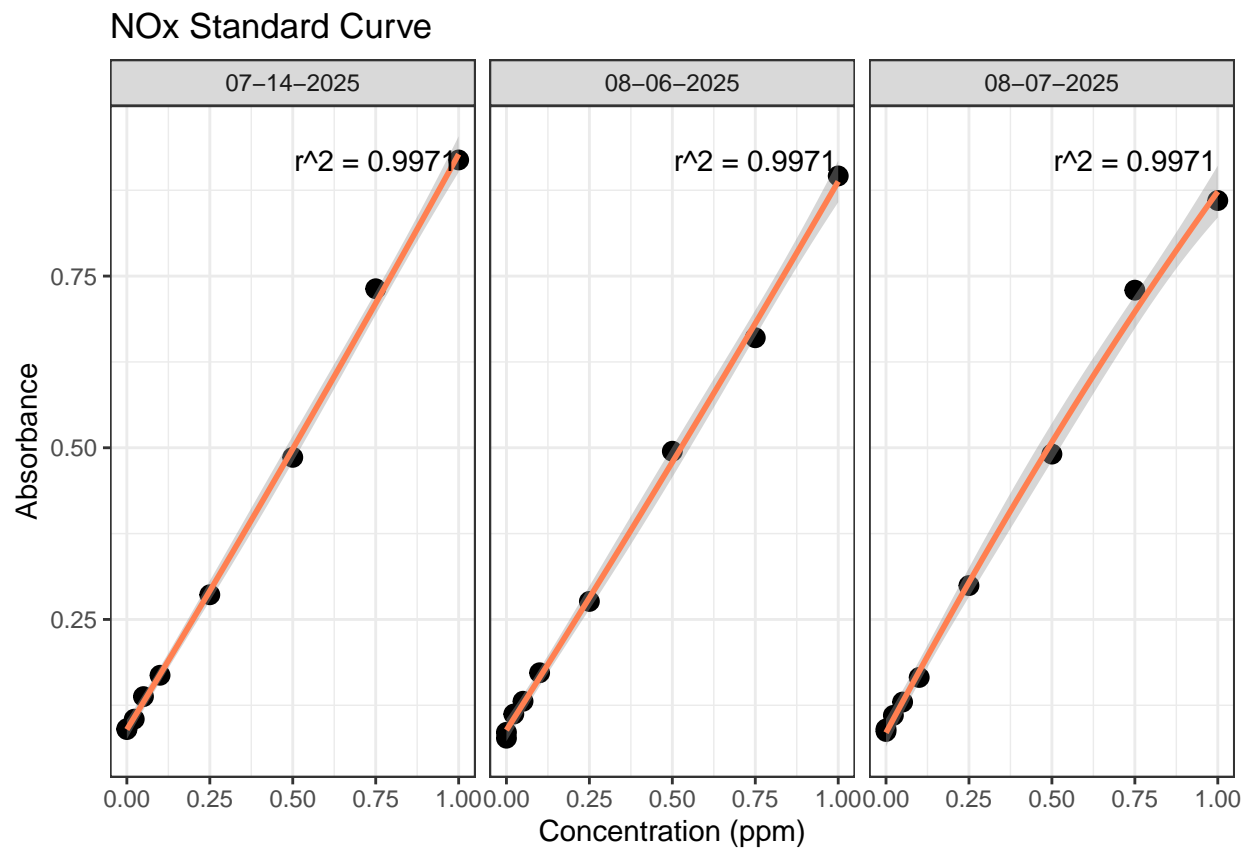


```
## 'geom_smooth()' using formula = 'y ~ x'
```

NH3 Standard Curve

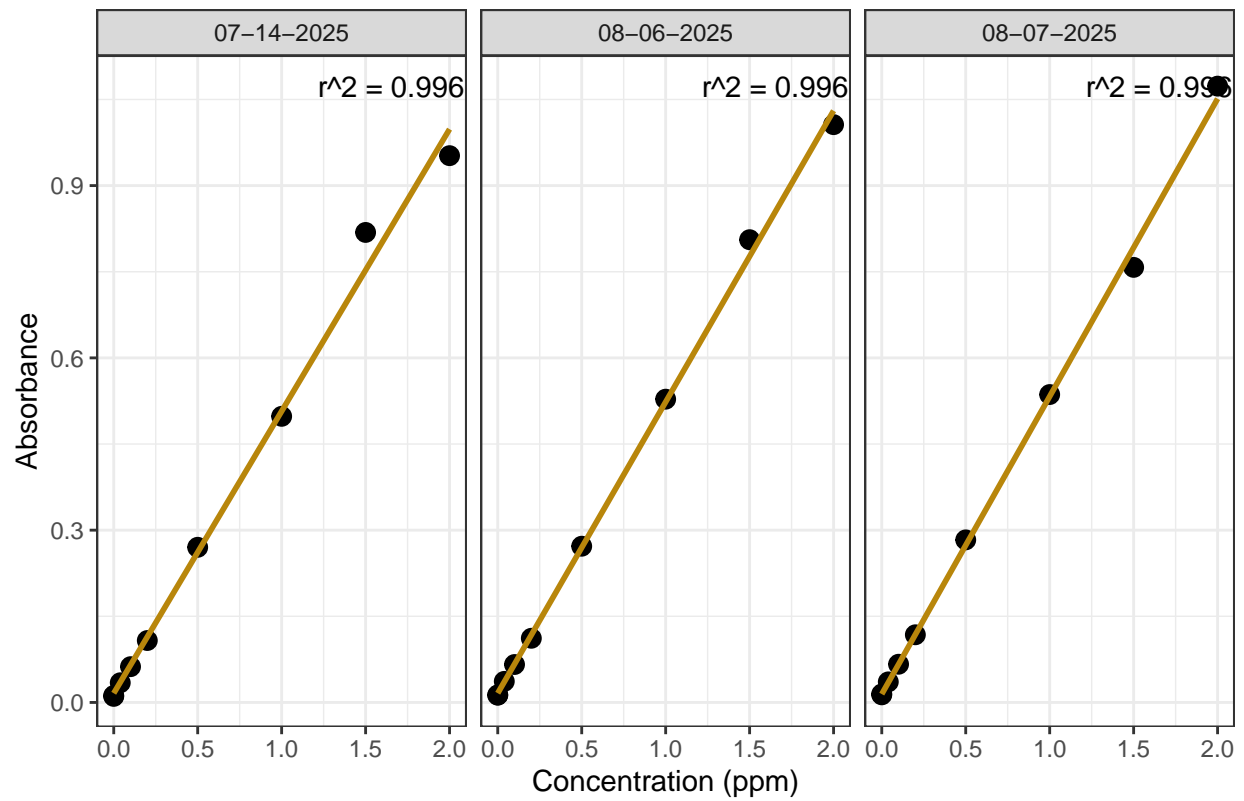

```
## 'geom_smooth()' using formula = 'y ~ x'
```

PO4 Standard Curve

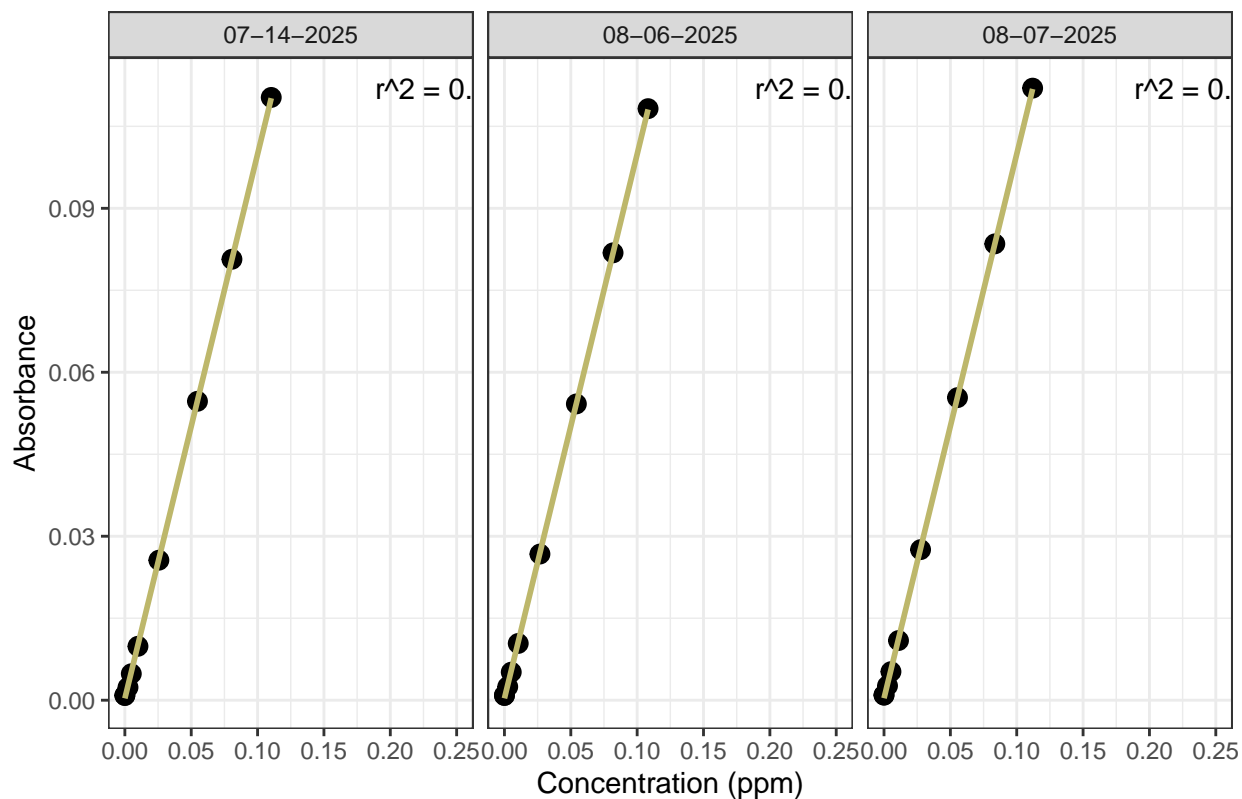

```
## [1] "NOx Curve r2 GOOD - PROCEED"
```

```
## [1] "NH3 Curve r2 GOOD - PROCEED"
```

```
## [1] "PO4 Curve r2 GOOD - PROCEED"
```

```
## [1] "QAQC log file exists and has been read into the code."
```

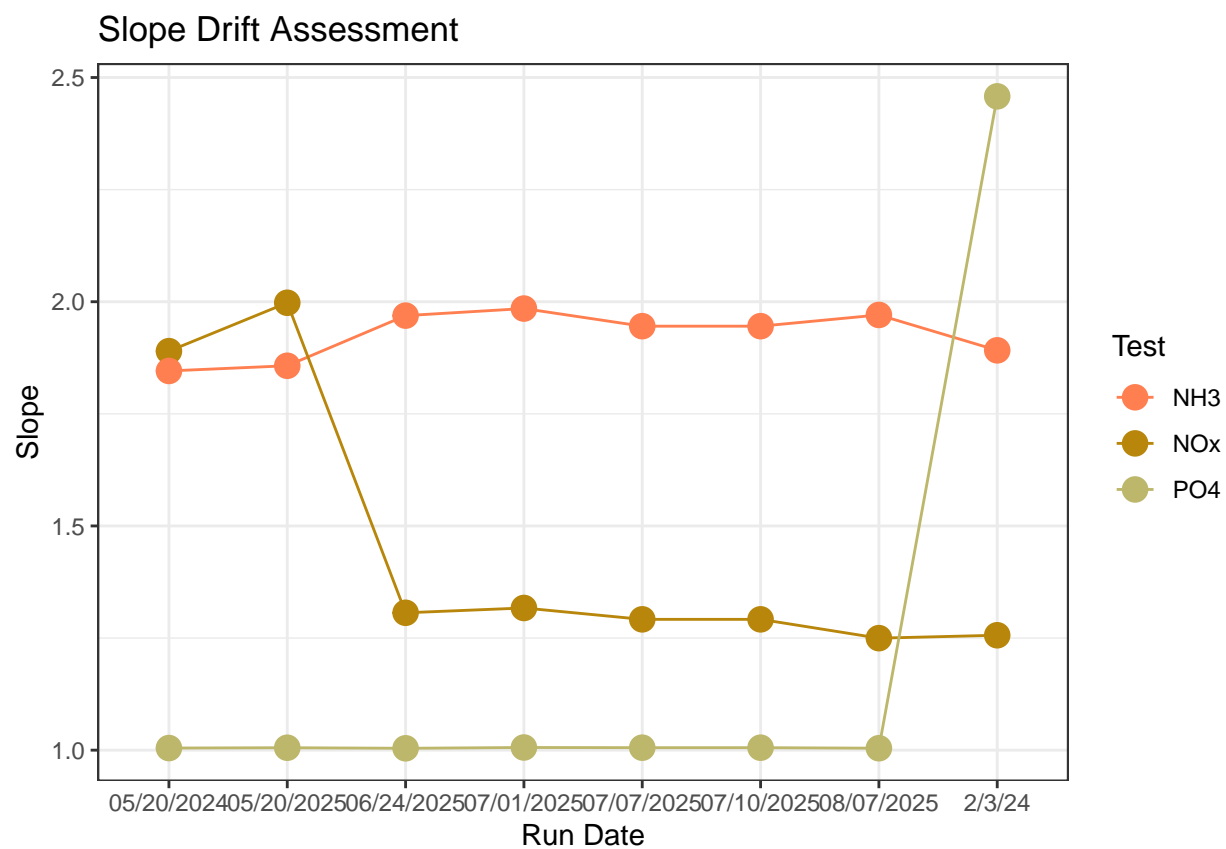


Table 1: Average Slope by Analyte

Test	avg_slope
NH3	1.926
NOx	1.450
PO4	1.187

0.3 Dilution Corrections - ensure the latest dilution is kept

```
## Dilution Corrections
```

```
## Duplicated samples: AbsChk_10ppt, AbsChk_20ppt, GCW_202411_WC_LysB_10cm, MSM_202411_TR_LysA_20cm, MS
```

```
##
```

```
## WARNING: Duplicated samples with NO dilution (possible input error):
```

```
## # A tibble: 6 x 4
```

```
##   Test      Sample_Name  n_dups any_dilution
##   <chr>      <chr>        <int> <lgl>
## 1 Ammonia 2    AbsChk_10ppt      2 FALSE
## 2 Ammonia 2    AbsChk_20ppt      2 FALSE
## 3 Vanadium NOx AbsChk_10ppt      2 FALSE
## 4 Vanadium NOx AbsChk_20ppt      2 FALSE
## 5 o-PHOS 0.3   AbsChk_10ppt      3 FALSE
## 6 o-PHOS 0.3   AbsChk_20ppt      3 FALSE
```

0.4 Performance Check

```
## [1] "NOx pe Check has a % Difference <25% - PROCEED"
```

```
## Run mean = 1.57518
```

```
## Expected  = 1.51
```

```
## [1] "NH3 pe Check has a % Difference <25% - PROCEED"
```

```
## Run mean = 1.209981
```

```
## Expected  = 1.034
```


```
## [1] "PO4 pe Check has a % Difference <25% - PROCEED"
```

```
## Run mean = 0.8969137
```

```
## Expected  = 0.824
```

```
##Check NOx Reduction Efficiency
```

```
## Assess Reduction Efficiency
```

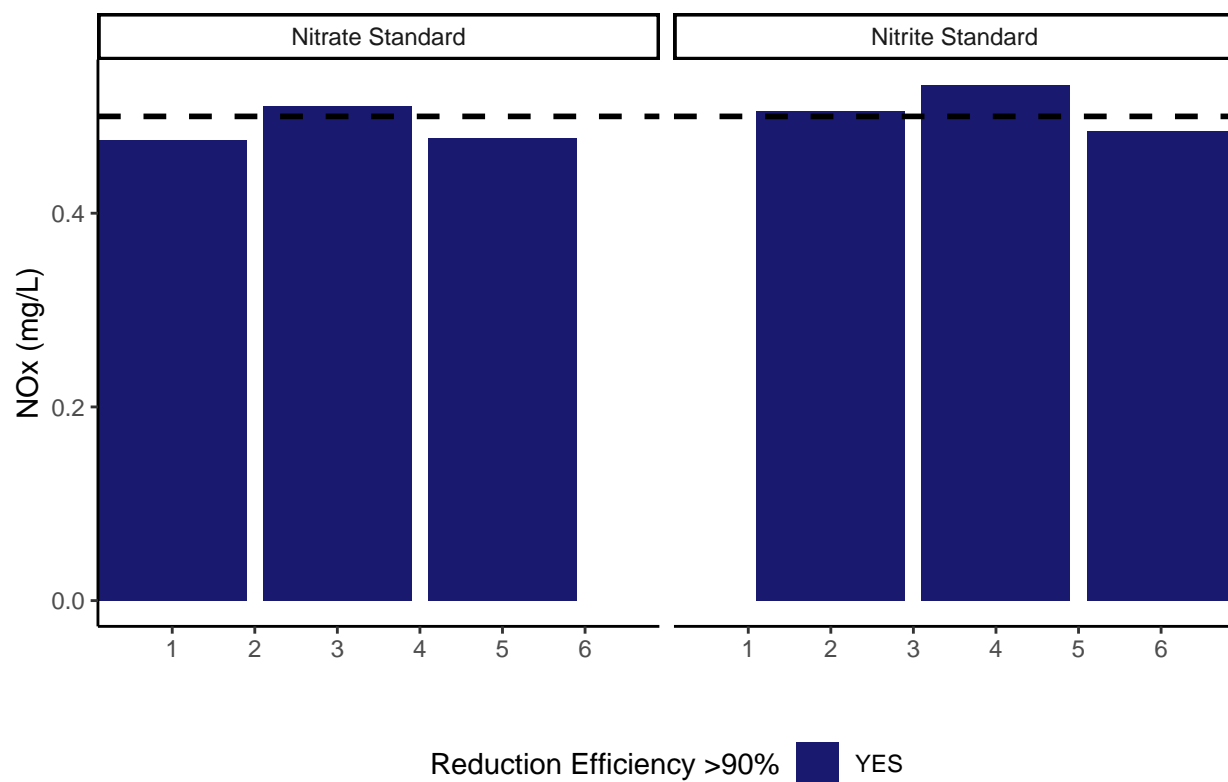

```
## [1] "Mean NOx Reduction Efficiency >95% - PROCEED"
```

```
## [1] 99.48873
```



0.5 Analyze the Check Standards

```
## Analyze Check Standards
```

```
## [1] "NOx Check Standard RSD within Range - PROCEED"
```

```
## [1] "NH3 Check Standard RSD within Range - PROCEED"
```

```
## [1] "PO4 CHECK STANDARD RSD TOO HIGH - REASSESS"
```

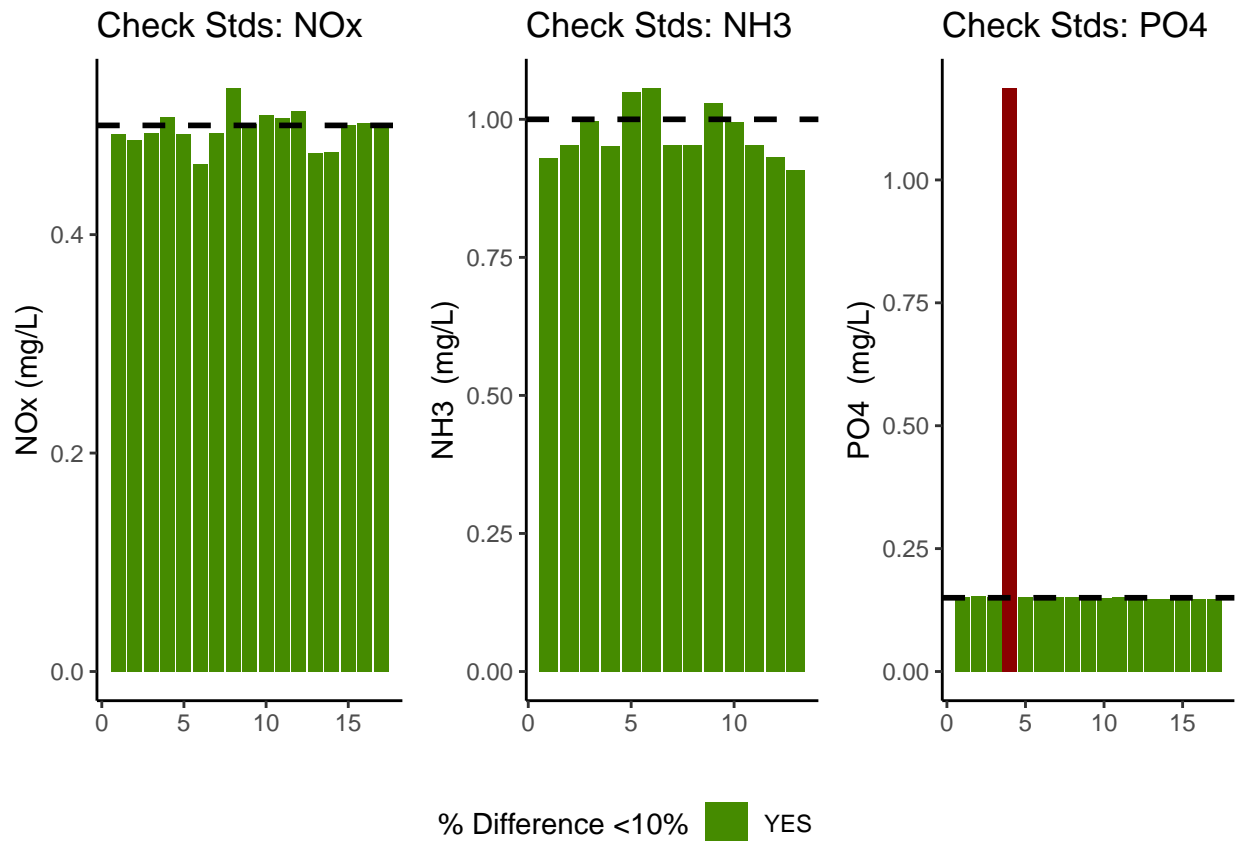

```
## [1] ">60% of NOx Check Standards are within range of expected concentration - PROCEED"
```

```
## [1] ">60% of NH3 Check Standards are within range of expected concentration - PROCEED"
```

```
## [1] ">60% of PO4 Check Standards are within range of expected concentration - PROCEED"
```

0.6 Analyze Blanks

```
## Assess Blanks
```

```
## [1] ">60% of NOx Blanks are below the lower 25% quartile of samples or 1/2 detection limit - PROCEED"
```

```
## [1] ">60% of NH3 Blanks are below the lower 25% quartile of samples - PROCEED"
```

```
## [1] ">60% of PO4 Blanks are below the lower 25% quartile of samples- PROCEED"
```

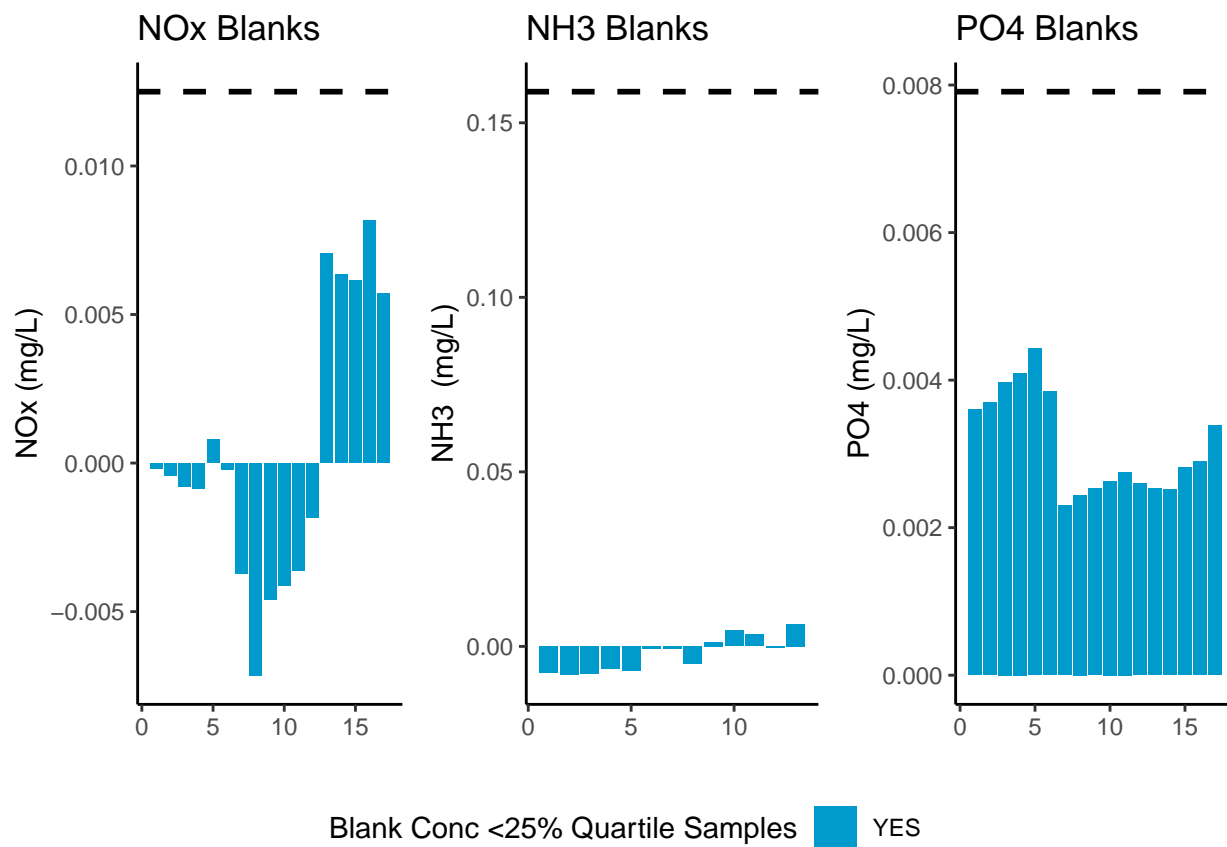


Table 2: Mean Concentration of Blanks

Test	Blank_Mean_Conc
NOx	0.0004
NH3	-0.0022
PO4	0.0031

0.7 Analyze Duplicates

```
## Analyze Duplicates
```

```
## [1] ">60% of NOx Duplicates have a CV <10% - PROCEED"
```

```
## [1] ">60% of NH3 Duplicates have a CV <10% - PROCEED"
```

```
## [1] ">60% of PO4 Duplicates have a CV <10% - PROCEED"
```


```
## Warning: Using 'size' aesthetic for lines was deprecated in ggplot2 3.4.0.
```

```
## i Please use 'linewidth' instead.
```

```
## This warning is displayed once every 8 hours.
```

```
## Call 'lifecycle::last_lifecycle_warnings()' to see where this warning was
```

```
## generated.
```

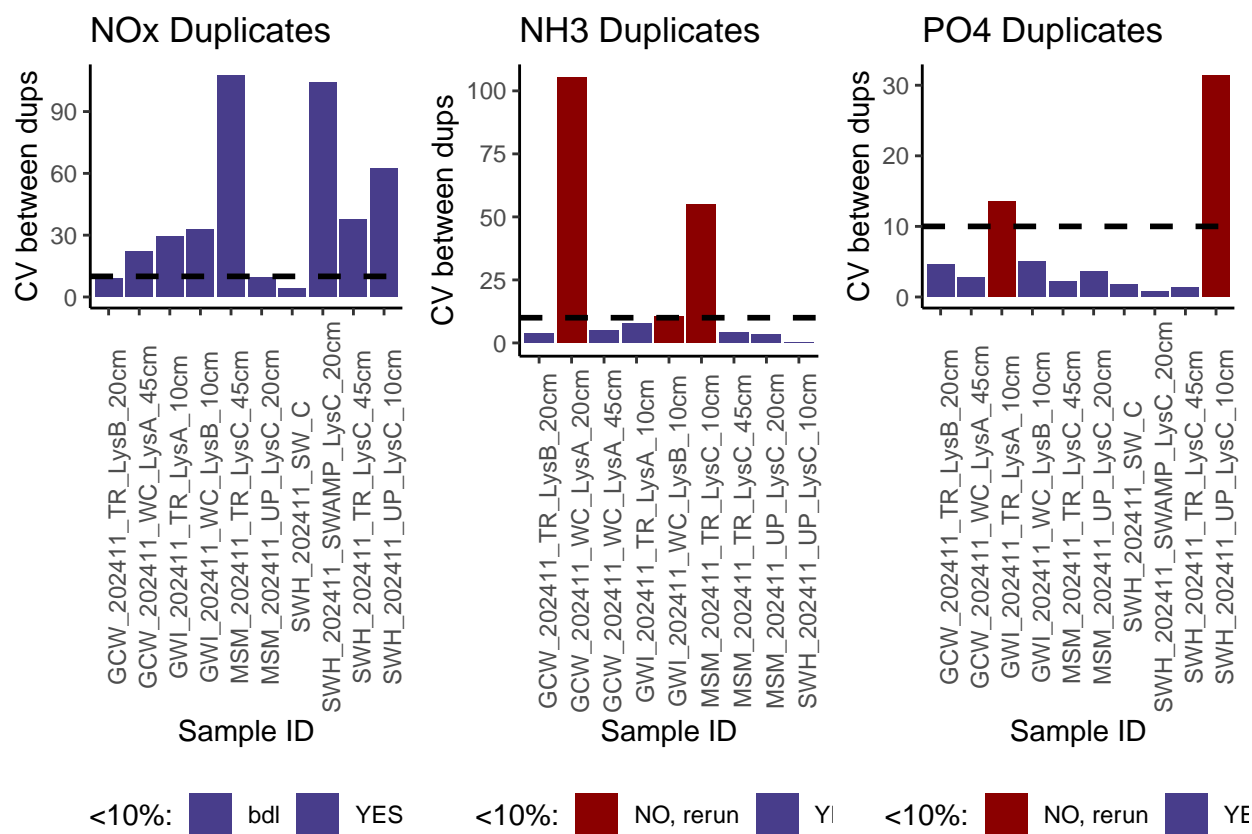


0.8 Spikes

[1] ">60% of Spikes have a CV <50% - PROCEED"

[1] ">60% of Spikes have a CV <50% - PROCEED"

[1] ">60% of Spikes have a CV <50% - PROCEED"

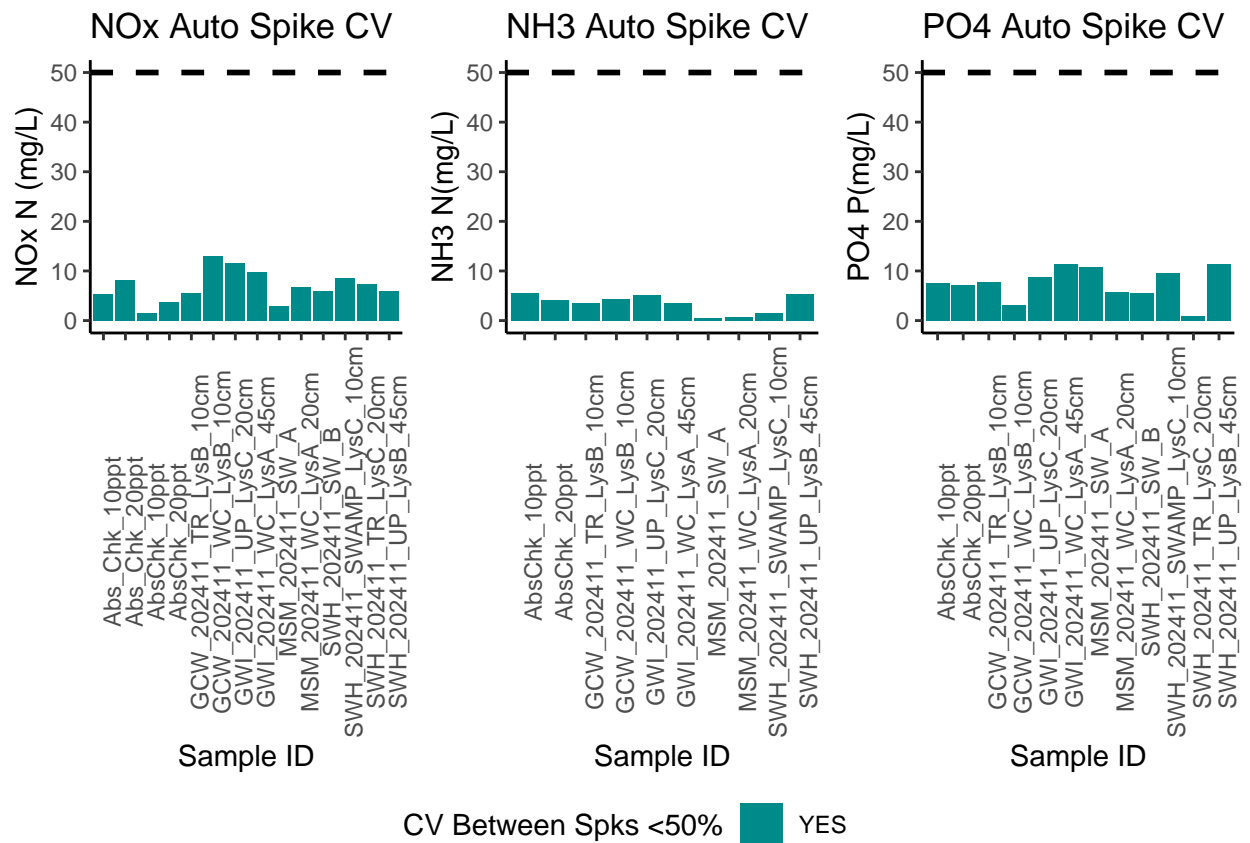

0.9 Matrix Effects

```
## [1] "NO NOx Matrix Effect, PROCEED"
```

```
## [1] "NO NH3 Matrix Effect, PROCEED"
```

```
## [1] "NO PO4 Matrix Effect, PROCEED"
```

0.10 Unit Converted Data Column Added (mg/L to uM)

0.11 Sample Flagging - Within range of standard curve

```
## Sample Flagging
```

0.12 Pull out sample id information

```
## Sample Processing
```

```
## Warning: Expected 5 pieces. Missing pieces filled with 'NA' in 33 rows [35, 36, 37, 44,  
## 59, 60, 70, 71, 72, 93, 94, 95, 137, 138, 139, 163, 178, 179, 188, 189, ...].
```

0.13 Check to see if samples run match metadata & merge info

```
## Check Sample IDs with Metadata
```


```
## Some sample IDs are missing from metadata.
```

```
## [1] "SWH_202411_SWAMP_LysC_45cm" "SWH_202411_TR_LysC_10cm"
```


```
## [3] "MSM_202411_TR_LysC_10cm"      "MSM_202411_TR_LysC_45cm"
```


```
## Warning: Expected 5 pieces. Missing pieces filled with 'NA' in 33 rows [35, 36, 37, 44,  
## 59, 60, 70, 71, 72, 93, 94, 95, 137, 138, 139, 163, 178, 179, 188, 189, ...].
```

Visualize Data

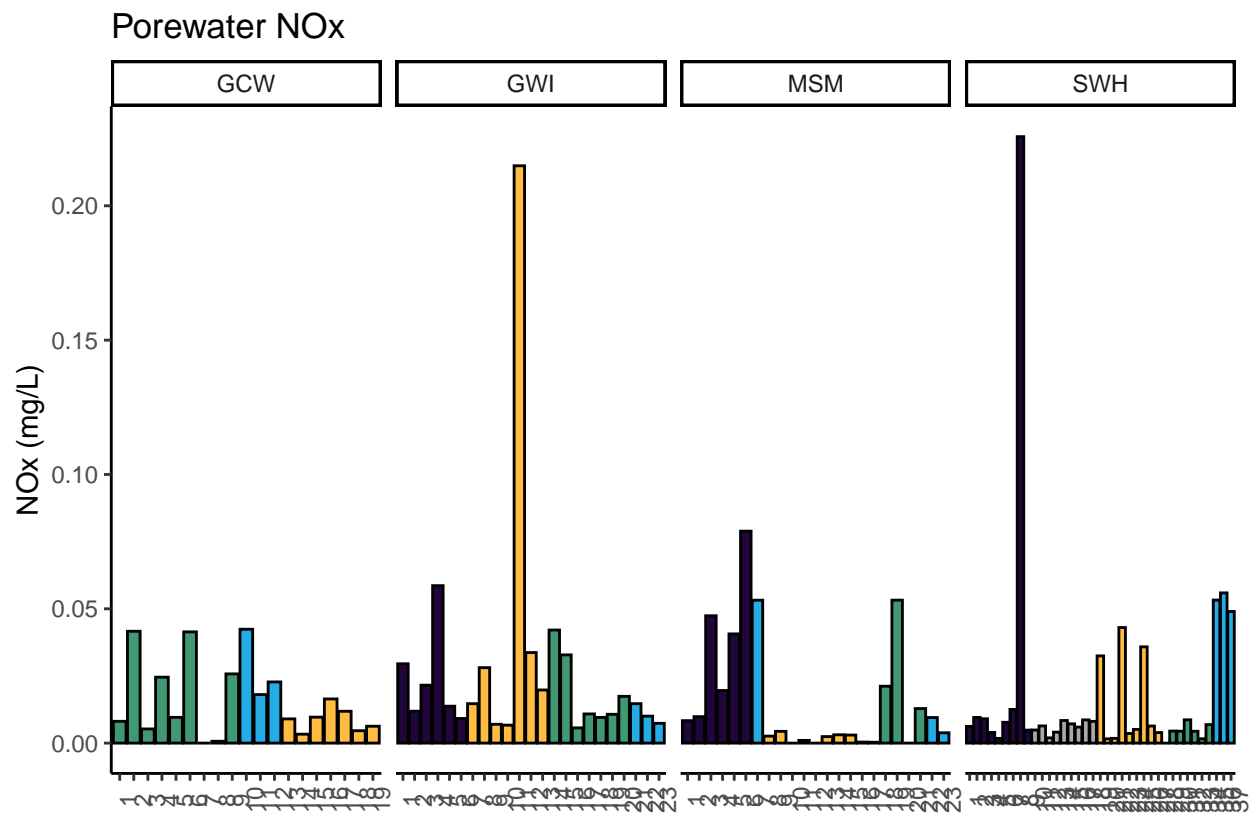

Porewater NH3

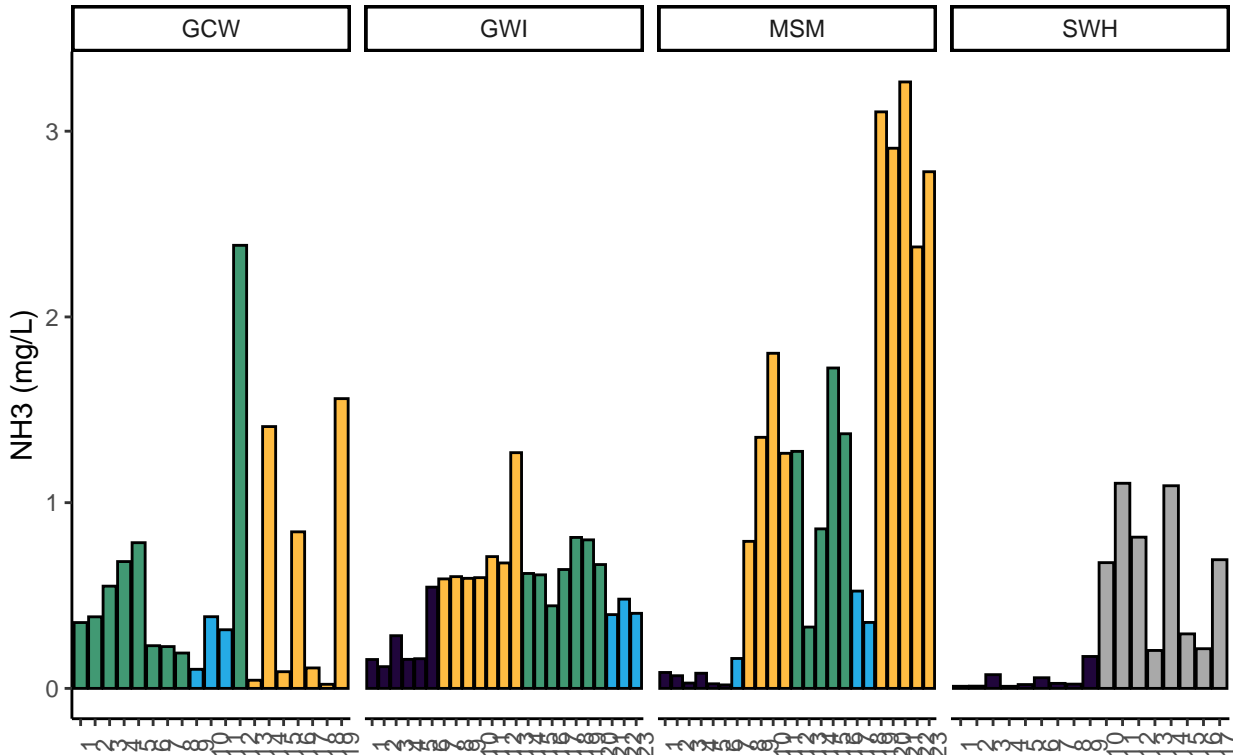

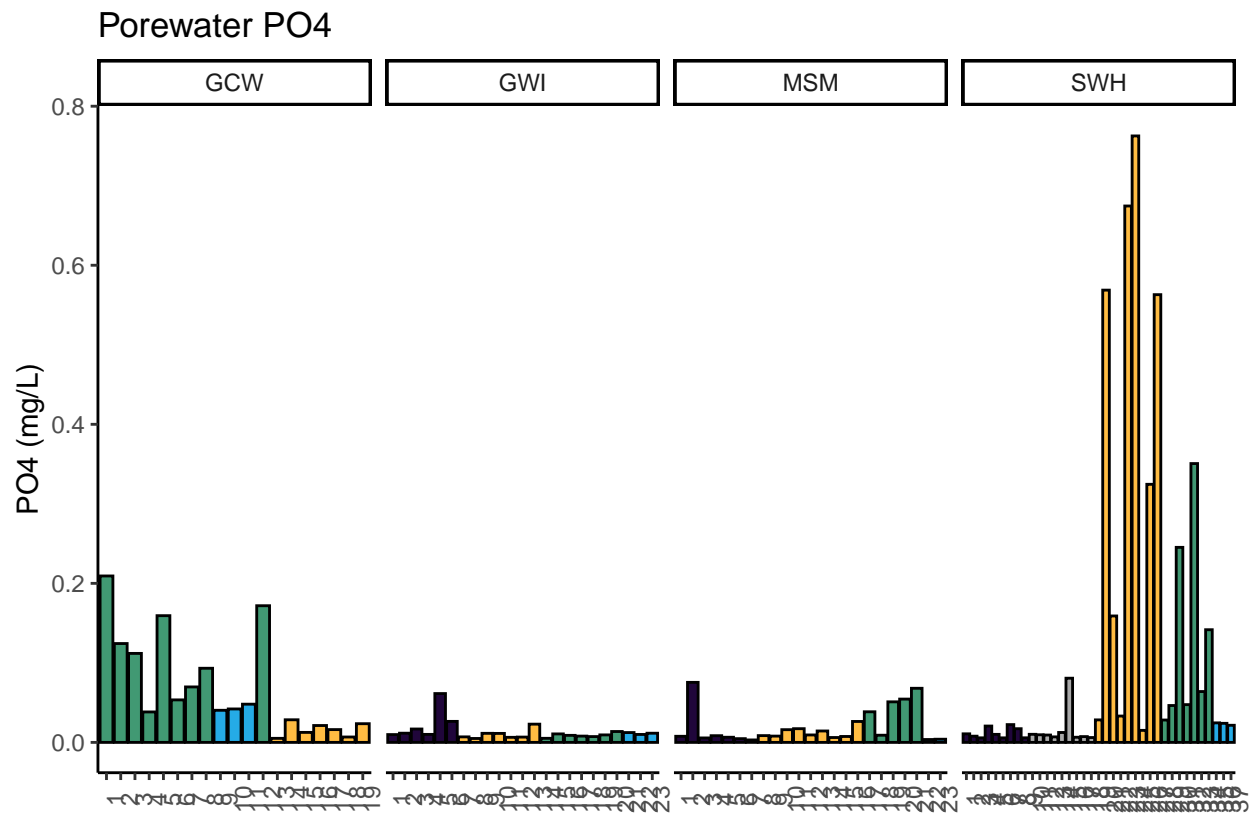

0.15 Export Processed Data

#end